




PP Fume Hood



LED Display

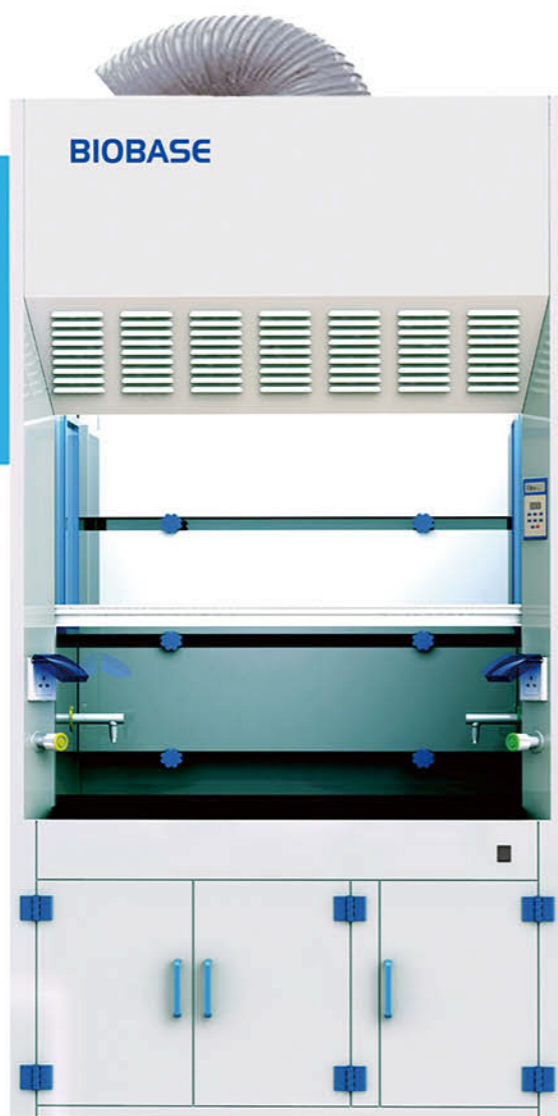

Waterproof Socket



PP Sink
resistant to acid,
alkali and anti-corrosion



Water&Gas
Remote Control

ADVANTAGE

1. It is safer to use anti-corrosive water tap.
2. Microprocessor control system, LED display.
3. With memory function in case of power-failure.
4. Made of porcelain white PP, resistant to acid, alkali and anti-corrosion.
5. Front window which is made of thick transparent toughened glass maximize light and visibility inside the fume hood, providing a bright and open working environment.

TECHNICAL PARAMETERS

Model	FH1000(P)	FH1200(P)	FH1500(P)	FH1800(P)
External Size (W*D*H)	1040*750*2200mm	1240*800*2200mm	1540*800*2200mm	1840*800*2200 mm
Internal Size (W*D*H)	820*520*872mm	1020*570*872mm	1320*570*872mm	1620*570*872 mm
Work Surface Height	740mm			
Max Opening	815mm			
Air Velocity	0.4~0.6m/s			
Noise	≤60dB			
Illuminating Lamp	LED lamp	Fluorescent Lamp	Fluorescent Lamp	Fluorescent Lamp
	12W*1	30W*1	30W*2	36W*2
Blower	Built-in PP axial flow blower (2 blowers for FH1800(P) only) ; Speed adjustable.			
Front Window	Resistant to acid, alkali Manual, 5mm toughened glass, height adjustable.			
Power Supply	AC220V±10%, 50/60Hz; 110V±10%, 60Hz			
Consumption	330W	360W	360W	360W
Material	Main Body	Made of Porcelain white PP, thickness 8mm, resistant to strong acid, alkali and anti-corrosion		
	Work Table	Chemical resistant phenolic resin		
Standard Accessory	Illuminating lamp, Water tap, Gas tap, Water Sink, Base cabinet Waterproof socket*2, 4 meters PVC duct (2 Pcs of 4 meters PVC duct for FH1800(P) only), Diameter: 250mm, PP centrifugal blower pipe strap*2 (4 Pcs for FH1800(P) only)			
Optional Accessory	PP work table, epoxy resin board or ceramic board Outer PVC centrifugal blower (only need to select when the duct more than 4m)			
Gross Weight	160kg	198kg	225kg	259kg
Package Size (W*D*H)	Main Body	1120*840*1620mm	1320*880*1620mm	1620*880*1620mm
	Base Cabinet	1120*840*1040mm	1320*880*1040mm	1620*880*1040mm