


MicroPad Shigella Count Plate

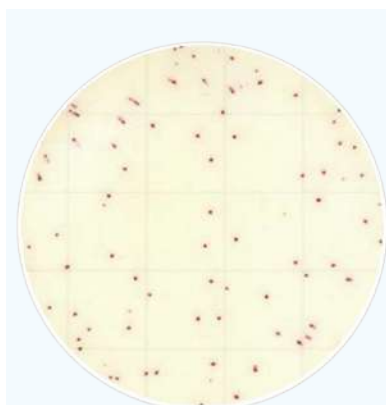


Product Code: MCPG02525 25 plates / pack

Shigella Count Plate (SHCP) is a sample-ready-culture medium system which contains readymade dry medium, a cold-water-soluble gelling agent and indicator. It can be used in direct counting of Shigella colonies after 24-48h incubation. The result is consistent with the corresponding ISO standards and other commercial counting plates.


**Read the user manual carefully before test, and follow the instructions.
Failure to do so may lead to inaccurate results.**




Sample Preparation

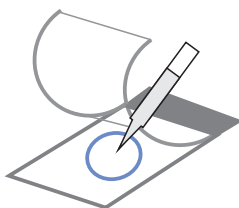

Sampling

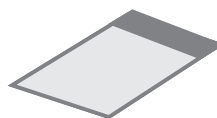

Incubation

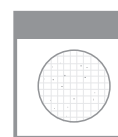

Enumeration

This Shigella Count Plate (SHCP) can be used in the quantitative plating of Salmonella Colonies in the food and beverage industries. Other applications are available upon request. Contact your supplier for more details.

Quantitative counting of the plate can be done by a standard colony counter or by software. All red colonies shall be counted regardless of the size or intensity.

Count all the red colonies.