


## MicroPad Yeast & Mould Count Plate

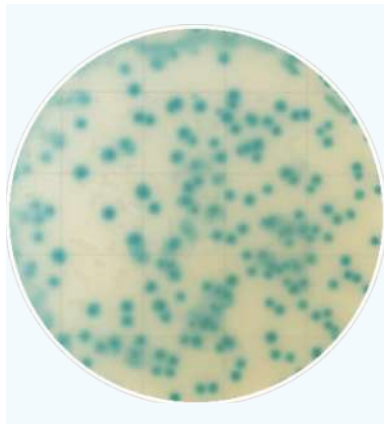


**Product Code: MCPG00825    25 plates / pack**

Yeast Mould Count Plate (YM) is a sample-ready-culture medium system which contains ready-made dry medium, somatomedin, a cold-water-soluble gelling agent and indicator. It can be used in direct counting of yeast and mould colonies after 2-3 days' incubation. The result is consistent with the corresponding ISO standards and other commercial counting plates.


**Read the user manual carefully before test, and follow the instructions.  
Failure to do so may lead to inaccurate results.**




**Sample Preparation**

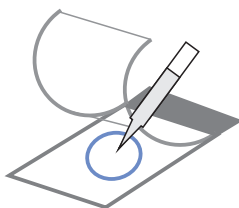

**Sampling**

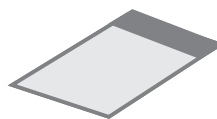

**Incubation**

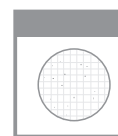

**Enumeration**

This Yeast Mould Count Plate can be used in the quantitative plating of yeast and mould colonies in the heat-processed product.

Quantitative counting of the plate can be done by a standard colony counter or by software. All green colonies shall be counted regardless of the size or intensity.

Count colonies within 10-150 CFU, count all the green colonies.