
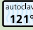





PharmaScoop



This sample scoop PharmaScoop has been developed in order to meet the high requirements of the pharmaceutical industry, such as those contained in the FDA regulations. Like all samplers for the pharmaceutical industry, they are made of stainless steel AISI 316 without edges, crevices or undercuts in which sample residues could collect and contaminate subsequent samples. The highly polished inner and outer surfaces of the scoop therefore meet all the requirements for GMP compliance in production and sampling.

- ▶ Stainless steel AISI 316 (1.4404)
- ▶ GMP conform
- ▶ Highly polished

PharmaScoop						
Capacity ml	Dia. mm	Length mm	Item no.			
10	25	45	5324-0010 N4			
50	38	70	5324-0011 N4			
100	50	100	5324-0012 N4			
500	95	140	5324-0018 N4			
750	95	210	5324-0013 N4			
1000	110	235	5324-0014 N4			
1250	125	195	5324-0015 N4			
1700	135	230	5324-0016 N4			
3000	185	215	5324-0017 N4			
PharmaScoop with long handle						
Capacity ml	Dia. mm	Length mm	Total length mm	Item no.		
1000	110	235	835	5324-0019 N4		





EasyScoop



When sampling bulk goods with a high specific weight, conventional scoops subject the user's wrist to great stress. EasyScoop with its centrally arranged inner grip is ergonomically designed and ensures easy and safe work. As EasyScoop is considerably shorter than conventional scoops, it is easier to take samples from narrow containers.

- ▶ Made of light aluminium, bright-anodised

EasyScoop					
Capacity ml	Length mm	Weight g	Item no.		
1500	200	325	5356-1500 N2		
2000	300	575	5356-2000 N2		